

| Service Note | | |
|--------------|---------------------------------|------------------|
| Product | Description: | Service note no. |
| BH250 | Spacers not fastened sufficient | SNo-T0174 |

Reason for note:

We have for a couple of production batches of BH250 sold in last half-year 2012 discovered 2 spacers and 2 screws which may not be fastened sufficient. This may cause a safety issue.

These production batches are within the following serial number ranges:


12403465 to 12405799

12404375 to 12404895 (Sold in the US)

BUT please notice that products in these ranges sold after January the 21st, 2013 are checked and OK.


Actions:

Make sure that the 2 spacers are securely fastened (torque = 65 Ncm)

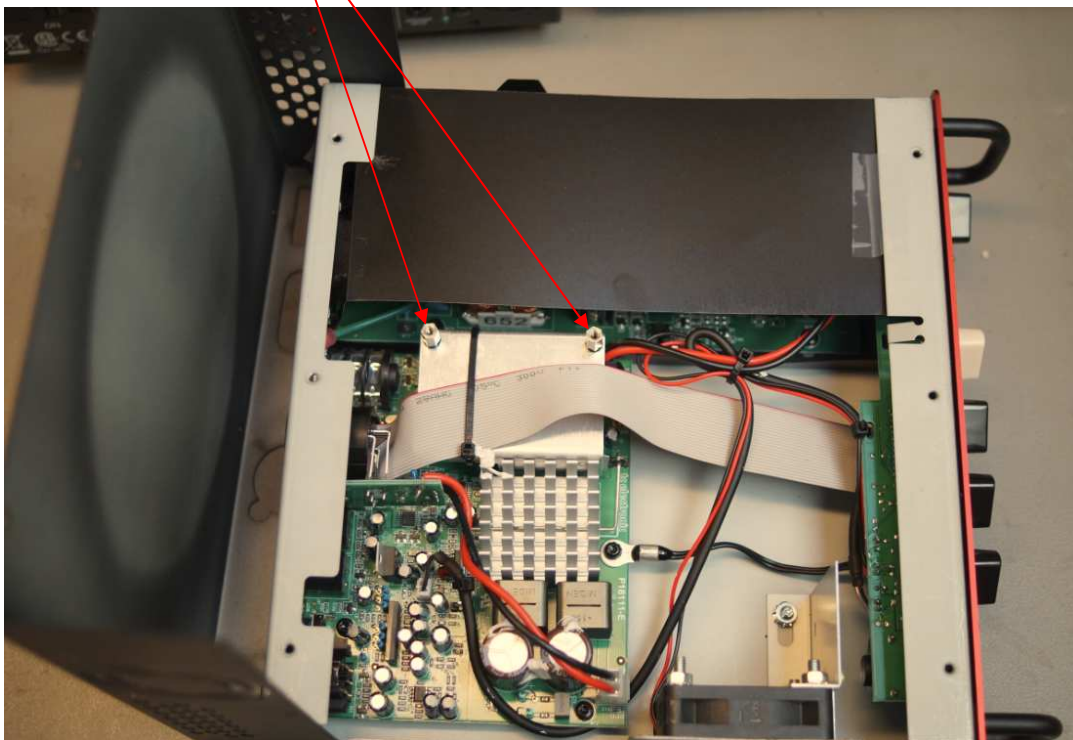

| Service Note | | |
|---------------------|--|------------------|
| Product | Description: | Service note no. |
| BH250 | Spacers not fastened sufficient | SNo-T0174 |

To access the 2 screws necessary to check, you have to take out the main board as shown. Make sure that the 2 screws are securely fastened (torque = 65 Ncm)

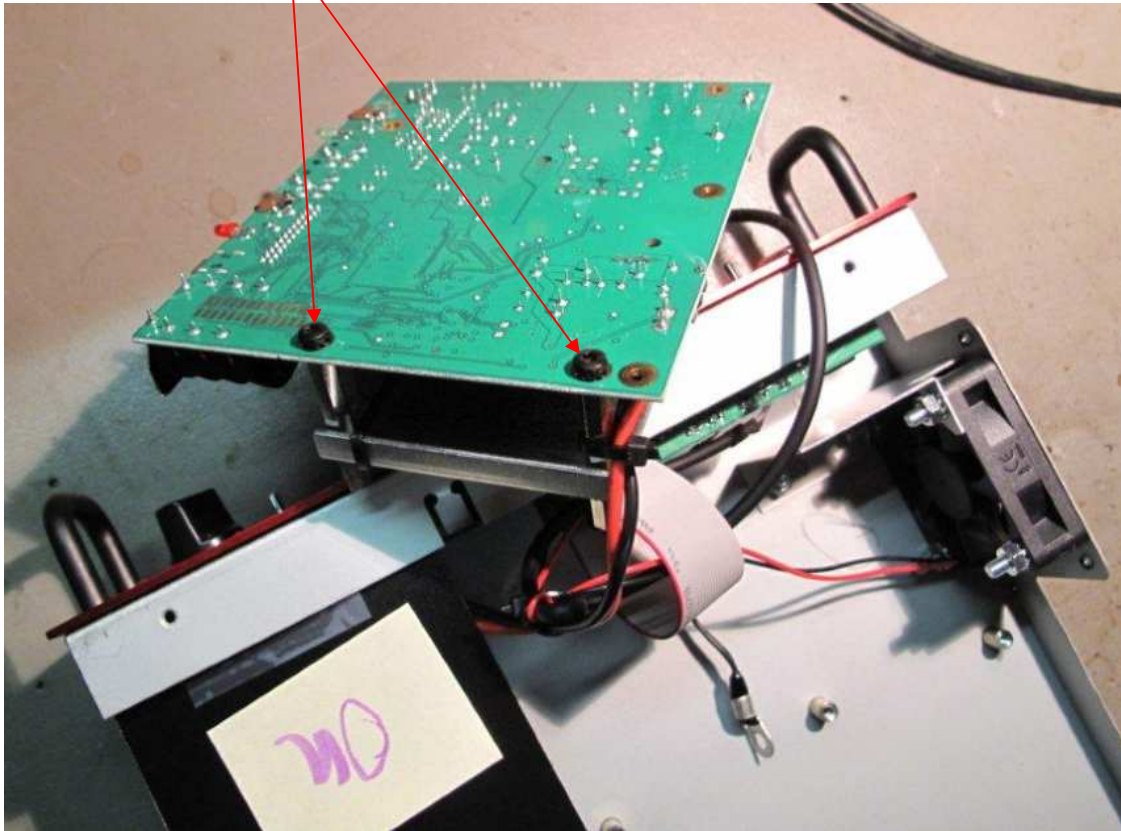

After making sure the screws and spacers are fastened correct, the main board shall be mounted again. Please pay attention that the mylar foil is placed correct.

It is recommended to place a paint mark dot or similar on the aluminium plate close to the spacers to indicate that this unit has been checked.

IMPORTANT:

This repair should only be carried out by skilled technicians.